


dataLink™ pro

Petri dish labeling solution

Ref. 439 050

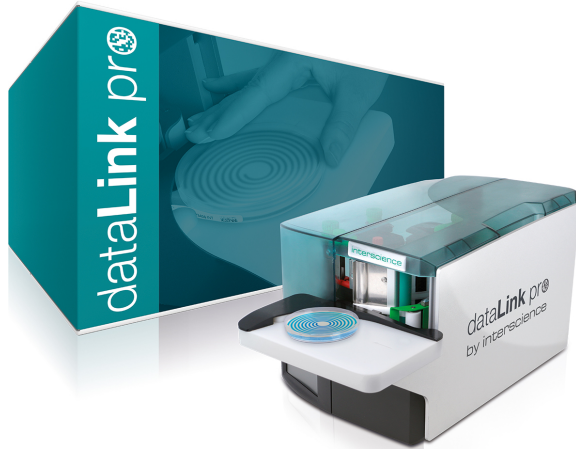

in compliance with:



Traceability:



Manufacture :



The dataLink™ pro pack is a semi-automated labeling solution for Petri dishes. Thanks to the dataLink™ software, the user manages and prints the information of his choice on a label which is stuck onto a Petri dish.

- After incubation, this data, from an easySpiral® plater, configures the automatic counting of colonies on Petri dishes by a Scan® counter (volume, dilution, time, deposit mode, etc.).
- Imported from an external file, this data will generate a series of labels which will be printed according to a template designed by the operator.

DATALINK™ SUPERVISION AND PRINTING SOFTWARE

- User-friendly interface in English
- Complies with 21 CFR Part 11 regulations requiring the control of systems handling digital data, in order to guarantee their integrity and security.
- User account management (connections, authorizations, action history)
- View audit trail or global activity log
- Monitoring of plating parameters processed by easySpiral®, Scan® and ScanStation® devices
- Automatic label printing, triggered by easySpiral® plating
- Import of XLSX, CSV or TXT files with data recognition
- Import of raw data by barcode scanning
- Customized labeling according to a print template
- Share Excel™ files with LIMS or external database

DATALINK™ PRO PRINTER

- Thermal transfer printing
- Chuck diameter: 25 mm
- Print resolution: 600 dpi
- Label size: 50 x 5 mm, thickness: 0.03 - 0.6 mm
- Label detection cell with position indicator
- Centered material guidance
- Semi-automatic labeling tray
- Pre-peel labeling
- USB 2.0 Hi-Speed port: computer connection
- Front USB port
- LCD display

1 / 2

dataLink pro_FT 11/23 Information and photos are not contractual. INTERSCIENCE reserves the right to change or improve product specifications without notice. Please visit www.interscience.com for more information. RCS 950 356 220 Versailles. INTERSCIENCE SARL F78860

interscience in the world

PARIS

Phone: +33 1 34 62 62 61
info@interscience.com

FRANKFURT

Phone: +49 611 7238 7770
sales.germany@interscience.com

BOSTON

Phone: +1 781 937 0007
sales.usa@intersciencelab.com

SHANGHAI

电话: +86 178 2123 6642
sales.china@interscience.cn

SINGAPORE

Phone: +65 6977 7232
sales.asia@interscience.com

TOKYO

Phone: +81 3 6712 9715
sales.japan@interscience.com

- Supply voltage and frequency: 100-240V ~ 50-60Hz
- Power: standby < 10 W / nominal 100 W
- Dimensions (W x D x H): 25.2 x 28.8 x 46 cm, net weight: 12 kg
- Box (W x D x H): 70 x 35 x 40 cm, gross weight: 16 kg

- Warranty: 1 year, upon registration of warranty card
- Optimum life of thermal transfer film and one roll of self-adhesive labels: 1 year
- Optimum label storage temperature: between 15°C and 25°C
- Made in Germany
- In compliance with EU regulations and Good Laboratory Practices (GLP)

DELIVERED WITH

- 1 power cord
- 1 USB cable
- 2 user manuals
- 1 thermal transfer ribbon
- 1 chuck
- 3 rolls of 2,500 self-adhesive labels, 1 roll mounted on printer
- 1 USB key with dataLink™ software

COMPUTER REQUIREMENTS

- Processor: Intel i5, 2.8 GHz, quad Core (or higher)
- Operating system: Windows™ 10 or 11 (or higher)
- RAM: 4 GB
- Equipment: 4 free USB ports (USB hub compatible)
- Screen resolution: 1280 x 1024 pixels (or more)

OPTIONS / ACCESSORIES

- Ref. 439 051 **Printer labels (dataLink™ pro): 3 rolls of 2,500 labels for dataLink™ printer pro - Label size: 50 x 5 mm - Quantity: 3 x 2500**
- Ref. 439 052 **Printer film (dataLink™ pro) : Thermal transfer film for dataLink™ pro printer**
- Ref. 439 170 **Petri dish automatic identification module (barcode (1D) and Datamatrix (2D) reader)**
- Ref. 522 000 **Barcode reader**